
## **Mindfulness - Application Flow**

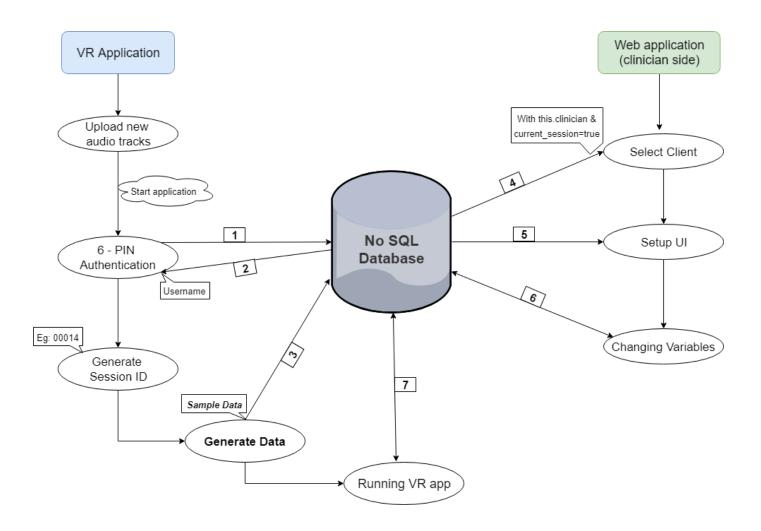

1	Once the client enters the 6-digit PIN, System authenticates with database
2	If authentication is successful, Database returns the username back to the VR app
3	The generated data is stored under the current username's specific session ID (created in the previous step), on the Database
4	Queries for clinician's clients and client's Current_session field (=true) and lists them on the web application
5	Queries the scene and audio settings from the generated data from (3) and displays in a neat UI
6	Changes made by the clinician on the scene & audio settings is stored into database OR retrieve changed variables from database, set from (7)
7	Fetches new settings from database, made by (6) and storing changed variables from VR app into database

## **Data Model (JSON format)**

```
<Username>: {
     "Current_session": "True/False",
     "Date": "",
     "Audio": {
         "Tracks": ["guide1.mp3", "guide2.mp3", "~NewlyAddedTracks"],
         "Current": "guide1.mp3",
         "Vocal volume": 0-100,
         "Background_volume": 0-100
     "Scene":{
         "Scene_id": 0,
         "Timeline": 0-24, //24 Hour format
         "Cloud_speed": 1-10,
         "Breathing": "True/False"
     },
     "Start_timestamp": "~DatTime format",
     "End_timestamp": "~DateTime format",
     "Duration": "~TimeSpan format",
     "Health monitor": {
         "Start_heartrate": 0,
         "End heartrate": 0,
         "Max_heartrate": 0,
         "Min_heartrate": 0,
         "Average_heartrate": 0,
         "File": <File_Upload_URL>
```

## **Sample Data (JSON format)**

```
'abhijit.93@hotmail.com":{
"00011":{
    "Current_session": "False",
"00012":{
    "Current_session": "False",
 "00014":{
    "Current session": "True",
    "Date": "",
    "Audio": {
        "Tracks": ["guide1.mp3", "guide2.mp3", "guide3.mp3", "myNewGuide.mp3"],
        "Current": "guide1.mp3",
        "Vocal volume": 70,
        "Background_volume": 50
    "Scene":{
        "Scene_id": 0,
        "Timeline": 0830,
        "Cloud_speed": 4,
        "Breathing": "True"
    "Start_timestamp": "3/12/2018 10:03:06 PM",
    "End timestamp": "NA",
    "Duration": "NA", //"02:14:18 TimeSpan format"
    "Health monitor": {
        "Start_heartrate": 0,
        "End_heartrate": 0,
        "Max_heartrate": 0,
        "Min_heartrate": 0,
        "Average_heartrate": 0,
        "File": "NA"
```